FitSpec

USER MANUAL FS-01

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FitSpec

USER MANUAL

- ◆ Please read the manual carefully and use the product correctly as instructed.
- ◆ After reading, keep the manual for future reference.
- ◆ If the product is transferred to others, ensure the manual is included to guarantee safe use.

Catalog

Safety Precautions	3
Product Information	4
Initial Use and Measurement	4
Starting Measurements	5
Maintenance and Cleaning	8
Product Certification Standards	8
Specifications.....	9
Product Warranty Information.....	10

Safety Precautions


- Please read the safety precautions before using this product.
- Follow the instructions below carefully.

 Warning

 Strict Compliance Required	This product is designed for home use and is only for health record-keeping, not medical purposes.
	If you feel unwell while using this product, stop immediately and consult a doctor before continuing.
	Not recommended for use by pregnant women.
	Do not use if you have implanted electronic medical devices, as electromagnetic interference may occur.
	Recommended age for use is 18 and above; for infants, only measure carried weight.
	Do not step on the analyzer with wet feet to avoid the risk of slipping.
	Place the product on a flat, stable surface. Do not step forcefully, and maintain a steady posture during measurement. Using it on uneven surfaces, such as carpets or cushions, may lead to inaccurate results.
	Avoid placing the product in direct sunlight or near heat sources to prevent damage to electronic components.
	Use only the specialized batteries designed for this product. Do not use these batteries for other devices, and do not replace them yourself.
	This is a precision instrument. Avoid dropping, shaking, or impacting it with external force.
	Store the body fat analyzer in a safe, dry, and ventilated place. Avoid using it in high-humidity areas, such as bathrooms.
	Avoid contact with sharp objects to prevent scratches and damage.

Product Information

- A. Detection Handle
- B. Handle electrodes
- C. Display Screen
- D. Foot Electrodes
- E. Type-C Charging Port
- F. Anti-Slip Foot Pads

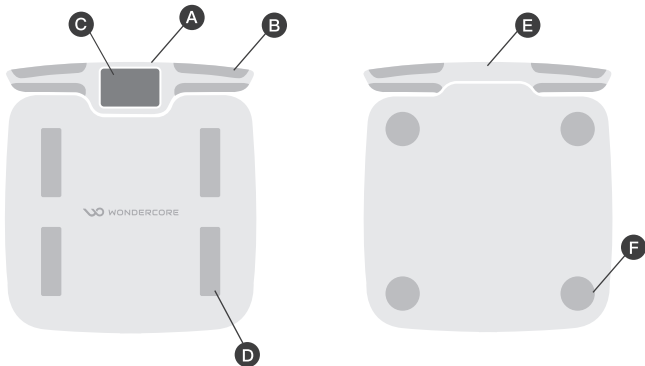Initial Use and Measurement


STEP 1 Download the Wonder Core App

1. Search for the app or scan the QR code to download it.
2. Follow the app instructions to register, log in, and complete your personal information settings.

STEP 2 Connect the Device

1. Place the scale on a flat surface.
2. Open the **Wonder Core App**, navigate to the device page. 
3. Tap the " + " symbol in the top right corner and select FitSpec.
4. Follow the app instructions to connect.

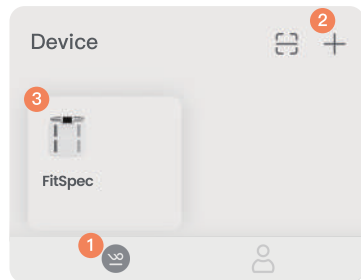STEP 3 Successful Binding and First Use

Once the device is successfully connected and the initial measurement is complete, tap "FitSpec" on the **device page**  to view results and historical records.

Starting Measurements

Method 1: Direct Measurement


1. Stand barefoot on the scale, placing both feet on the foot electrodes. The display screen will light up.
2. Pull up the detection handle, ensuring all fingers and palms are in contact with the handle electrodes.
3. Follow the screen prompts to begin measurement. Once completed, view results on the screen.

Note : This mode can only measure weight as it is not connected to the app.

Method 2: Measurement with App Connection

1. Log into the app, tap “FitSpec” , and select “Start Measurement” to wake the scale and connect.

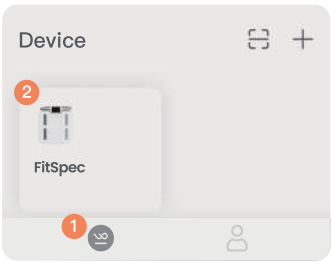
Note : Connection is only possible when the screen is illuminated.

2. Stand barefoot on the scale, placing both feet on the foot electrodes.
3. Pull up the detection handle, ensuring all fingers and palms are in contact with the handle electrodes.
4. Follow prompts to complete measurement. Results will be available in the app.

Reset Instructions

Press and hold the buttons on both sides of the handle simultaneously until the device restarts.

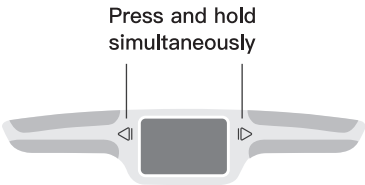
Posture During Measurement

To ensure the stability and accuracy of measurement results, maintain the following correct postures during measurement:


1. Torso and Limbs

- Keep your arms from touching your abdomen, arms should not be raised too high. Extend them straight, ensuring that the armpits do not touch the sides of your body.
- If the armpits touch the body, wear lightweight, sleeved clothing.

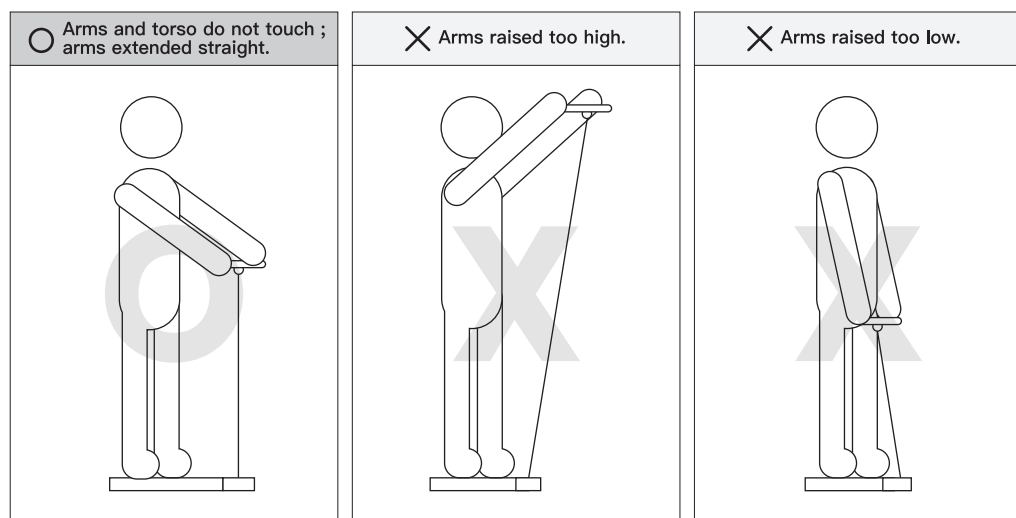
2. Hands

- Ensure both hands are separated and that fingers and palms are fully in contact with the handle electrodes.

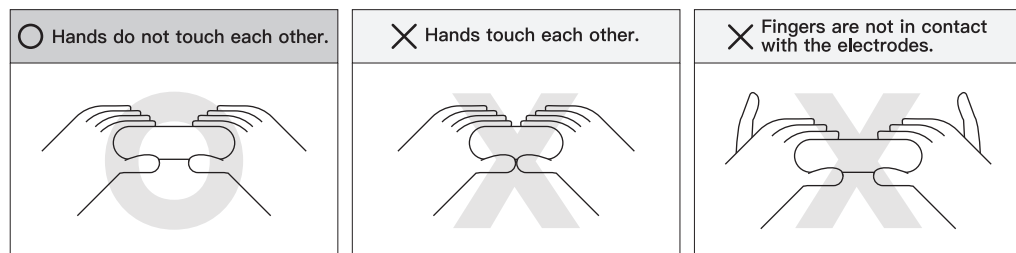3. Feet

- Keep both feet apart. If legs touch, wear lightweight shorts during testing.
- Ensure feet are fully in contact with the foot electrodes, barefoot. Socks are not allowed.

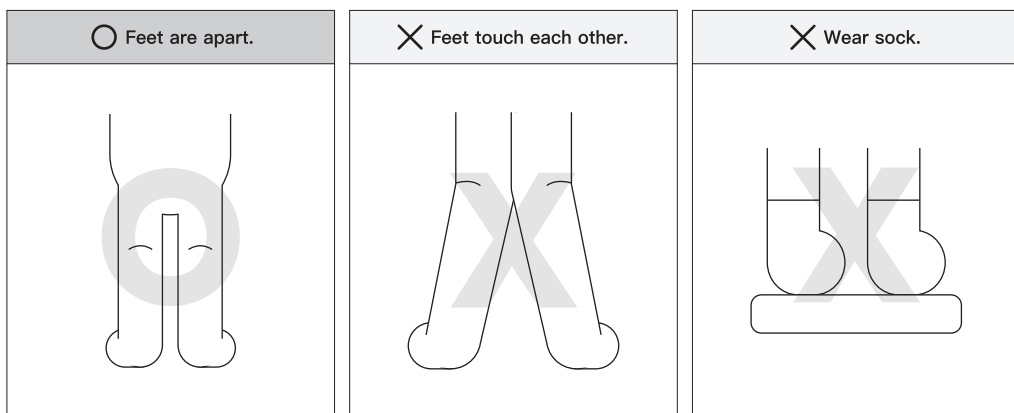Ensuring Accurate Measurements

To achieve accurate results, follow these guidelines to eliminate interference factors :

1 Input Accurate Personal Information :

Enter correct gender, age, and height.
Incorrect details can affect measurement results.

2 Avoid Post-Meal Measurements :

Food in the stomach may impact accuracy. Test on an empty stomach or wait 3~6 hours after eating.

3 Remove Metal Accessories :

Before measurement, remove metal jewelry and accessories.

4 Test in the Morning :

Long periods of standing or sitting can cause water retention in the lower body. Morning tests are recommended.

5 Avoid Post-Exercise Testing :

Exercise temporarily alters body composition, affecting results.

6 Avoid Testing After Bathing, Massages, or Soaking in Water :

Increased blood flow and sweating can alter body composition temporarily.

7 Stand for 5 Minutes Before Testing :

Prolonged lying or sitting affects water distribution in the body, impacting accuracy.

8 Wear Lightweight Clothing :

Avoid wearing heavy garments during measurement.

9 Minimize External Factors :

Conduct tests at a consistent time and location for reliable long-term statistics.

Maintenance and Cleaning

1. Regularly use a soft, dry cloth or alcohol wipes to gently clean the scale's surface, removing dust and dirt.
2. Avoid using strong acidic or alkaline cleaners, as they may damage the scale's surface or internal electronic components.
3. Prevent water or other liquids from entering the scale to avoid damage when cleaning.

Product Certification Standards

IC Compliance

Innovation, Science and Economic Development (ISED) Canada Regulatory Compliance

CAN ICES-3(B)/NMB-3(B)

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

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This equipment complies with Innovation, Science and Economic Development Canada's radiation exposure limits set forth for an uncontrolled environment. The End user must follow the specific operating instructions for satisfying RF exposure compliance. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

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FCC Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

RF exposure warning

The equipment complies with FCC RF exposure limits set forth for an uncontrolled environment. The equipment must not be co-located or operating in conjunction with any other antenna or transmitter.

Specifications

Model	WC-FS-01		
Dimensions (L x W x H)	375 x 352 x 35 mm		
Weight	2.7 kg		
Weight Limit	182 kg		
Measurement Range	0.1 ~ 180 kg	Measurement Units	kg / st:lb / lb
Battery Capacity	800 mAh	Charging Method	Type-C
Charging Time	2 ~ 2.5h	Usage Time	Approximately 90 days
Display Screen	3.5-inch TFT screen		
BLUETOOTH®	5.0		
Materials	ABS 、tempered glass 、stainless steel 、rubber		
Manufacturing Location	Patented in Taiwan, manufactured in China		
Purpose	Body fat scale		
Manufacturer/Importer	WONDER CORE CO., LTD.		
Address	32F, No. 213, Chaofu Rd., Xitun Dist., Taichung City 407, Taiwan		
Customer Service(Whatsapp)	1(954)243 2260		
E-mail	supportus@wondercore.com		
Website	www.wondercore.com		

 This product is for home use only and is intended for health record-keeping purposes, not for medical use.

Product Warranty Information

● Scope of Warranty :

Products manufactured, authorized, or sold by the company under legitimate sales and without payment disputes fall under this warranty's scope. To request warranty service, contact Wondercise Technology Corp. (the company) and present the original purchase invoice or a warranty card stamped by the company. Keep these documents safe for warranty claims.

● Warranty Content :

The manufacturer provides a one-year warranty, effective from the date of purchase. If the product has defects due to manufacturing or shipping, replacements are allowed within seven days of receipt. Warranty terminates if the product is resold or transferred to another owner.

During the valid warranty period, if functional defects occur under normal use due to manufacturing or assembly flaws in the product, the company will unconditionally repair or replace the defective parts free of charge. The original manufacturer has the right to inspect and handle the replaced parts under the warranty terms. All replaced parts will become the property of the company. The terms of this warranty also apply to any parts repaired or replaced free of charge during the warranty period.

● Warranty Process :

- (1) The customer identifies product defects or malfunctions.
- (2) Verify labels and packaging to confirm it is the company's product and within the warranty period.
- (3) Contact the company and describe the issue. Upon confirmation, the complaint is processed under warranty terms.
- (4) The company will arrange product recovery and compensation based on actual transactions.
- (5) After free repair or replacement, customers cannot request refunds, exchanges, or further compensation for the same issue.

● Technical Support for the Product :

Before the product is sold, the company provides detailed product descriptions, specifications, and usage or storage precautions. Customers are obligated to fully understand the product information provided by the company and use or store it correctly.

● Exclusions :

Under the following circumstances, our company may waive warranty liability :

- (1) The defective product provided by the customer is confirmed to not be from our company.
- (2) The defective product has been used beyond the warranty period.
- (3) The product has not been serviced or maintained by an authorized dealer of our company, or defects, malfunctions, or damage have occurred due to misuse, abuse, or unauthorized replacement or repair of parts.
- (4) The defective product has been confirmed not to have been serviced or maintained by the original manufacturer or an authorized dealer, or defects, malfunctions, or damage have resulted from misuse, abuse, or unauthorized replacement or repair of parts.
- (5) The defective product has an improper origin, such as being obtained illegally or through financial disputes.
- (6) Natural wear and tear of product components.
- (7) Malfunctions caused by force majeure events such as fire, flood, or earthquake.
- (8) The product is a personal item and has been unsealed or used.
- (9) Malfunctions caused intentionally by the customer or due to failure to use and store the product as instructed.

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